



A Publication
of Reliable Methods
for the Preparation
of Organic Compounds

Working with Hazardous Chemicals

The procedures in *Organic Syntheses* are intended for use only by persons with proper training in experimental organic chemistry. All hazardous materials should be handled using the standard procedures for work with chemicals described in references such as "Prudent Practices in the Laboratory" (The National Academies Press, Washington, D.C., 2011; the full text can be accessed free of charge at http://www.nap.edu/catalog.php?record_id=12654). All chemical waste should be disposed of in accordance with local regulations. For general guidelines for the management of chemical waste, see Chapter 8 of Prudent Practices.

In some articles in *Organic Syntheses*, chemical-specific hazards are highlighted in red "Caution Notes" within a procedure. It is important to recognize that the absence of a caution note does not imply that no significant hazards are associated with the chemicals involved in that procedure. Prior to performing a reaction, a thorough risk assessment should be carried out that includes a review of the potential hazards associated with each chemical and experimental operation on the scale that is planned for the procedure. Guidelines for carrying out a risk assessment and for analyzing the hazards associated with chemicals can be found in Chapter 4 of Prudent Practices.


The procedures described in *Organic Syntheses* are provided as published and are conducted at one's own risk. *Organic Syntheses, Inc.*, its Editors, and its Board of Directors do not warrant or guarantee the safety of individuals using these procedures and hereby disclaim any liability for any injuries or damages claimed to have resulted from or related in any way to the procedures herein.

These paragraphs were added in September 2014. The statements above do not supersede any specific hazard caution notes and safety instructions included in the procedure.

Organic Syntheses, Coll. Vol. 8, p.190 (1993); Vol. 67, p.48 (1989).

TELOMERIZATION OF ISOPRENE WITH DIALKYLAMINE: *N,N*-DIETHYLNERYLAMINE

[2,6-Octadien-1-amine, *N,N*-diethyl-3,7-dimethyl-, (*Z*)-]

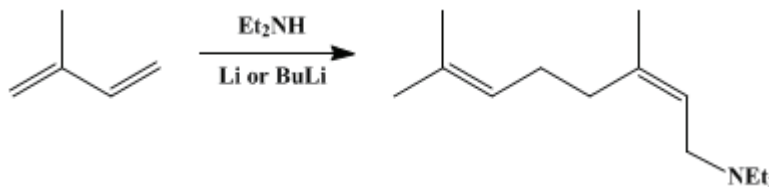

Submitted by Kunihiko Takabe¹, Takashi Yamada, Takao Katagiri¹, and Juntaro Tanaka¹.
Checked by Alan J. Chalk, Laszlo V. Wertheimer, and Gabriel Saucy.

1. Procedure

A 100-mL glass autoclave (Note 1) is charged with 34 g (0.50 mol) of isoprene (Note 2) and 7.3 g (0.10 mol) of diethylamine (Note 3) under nitrogen (Note 4). The contents are cooled to 3–5°C using an ice–salt bath, and 3.7 mL (0.006 mol) of a 1.62 *M* solution of butyllithium in hexane (Note 5) is added dropwise by syringe while stirring. The resulting slightly yellow solution is blanketed carefully with dry nitrogen, and the vessel is sealed. The solution is heated to 65°C and stirred for 13 hr. The vessel is then cooled using an ice–salt bath, the contents are poured into 15 mL of water, and the organic layer is separated. The aqueous phase is extracted with 30 mL of diethyl ether. The organic layers are combined, dried over anhydrous sodium sulfate, filtered, and concentrated under reduced pressure with a rotary evaporator. The residual liquid is distilled under reduced pressure (Note 6) to afford 16.1–18.0 g (77–86%) of *N,N*-diethylnerylamine as a colorless liquid, bp 83–84°C (1.5 mm) (Note 7) and (Note 8).

2. Notes

- The glass autoclave must be used with appropriate precaution. A pressure of about 20 psi is generated, so a low-positive-pressure reactor, such as a shielded Fischer-Porter Bottle (Fischer and Porter Company, Lab-Coest Division, Warminster, PA18974), can be used instead. Alternatively, the reaction may be carried out without pressure equipment in 3 days.²
- Isoprene, obtained from Aldrich Chemical Company, Inc. (also available from Nakarai Chemicals, Ltd. in Japan), was dried over Linde 4A molecular sieves for at least 1 day and freshly distilled prior to use.
- Diethylamine, obtained from Aldrich Chemical Company, Inc., was distilled from calcium hydride before use.
- A slow stream of dry nitrogen was passed through an inverted funnel that was placed over the vessel.
- Butyllithium in hexane was obtained from Aldrich Chemical Company, Inc. (also available from Nakarai Chemicals, Ltd. in Japan), and titrated before use.³ Lithium (0.055 g, 0.0079 g-atom), cut into small pieces (ca. 2 mm), can be used in place of the butyllithium-hexane solution.
- A forerun (1.4–0.6 g) boiling at 52–56°C (16 mm) consists mainly of *N,N*-diethyl-2-methyl-2-butenylamine and *N,N*-diethyl-3-methyl-2-butenylamine.
- GLC analysis (7% Apiezon L on 60–80-mesh Gaschrom Q, 3 mm × 3 m, 170°C) showed that the product is over 99% isomerically pure. The spectral properties of *N,N*-diethylnerylamine are as follows: IR (neat) cm⁻¹: 1660, 1200, 1165, 1050, 830; ¹H NMR (CDCl₃) δ: 1.03 (t, 6 H, *J* = 7, CH₃-), 1.45–1.7 (m, 6 H, CH₃-C=), 2.02–2.2 (m, 4 H, -CH₂-C=), 2.50 (q, 4 H, *J* = 7, -N-CH₂-), 3.05 (d, 2 H, *J* = 6.5, -N-CH₂-C=), 4.9–5.5 (m, 2 H, -CH=).
- The submitters report that similar reactions of isoprene with dimethylamine (55°C, 10 hr) and dipropylamine (65°C, 15 hr) afforded *N,N*-dimethylnerylamine (69%), bp 85–86°C (9 mm) and *N,N*-dipropylnerylamine (77%), bp 86–87°C (1 mm), respectively.

3. Discussion

The procedure described here is essentially that reported earlier⁴ and modified by subsequent experience.⁵ It is a simple, general method for the synthesis of an *N,N*-dialkylnerylamine from *isoprene* and a dialkylamine. *N,N*-Dialkylnerylamines can also be prepared by the reaction of a dialkylamine with *neryl chloride*⁶ by a modification of Sandler's method.⁷ However, pure *nerol*, the starting material of *neryl chloride*, is expensive, and not easily available commercially. The present method illustrates a mild and convenient one-step reaction for the preparation of an *N,N*-dialkylnerylamine. In addition, the starting materials are readily accessible, the reaction proceeds stereoselectively, and the yields of the product are generally high. This process consists of the initial addition of lithium dialkylamide to *isoprene*, followed by the propagation of the resulting intermediate, and the termination by dialkylamine.

N,N-Dialkylnerylamines serve as convenient precursors for naturally occurring acyclic monoterpenes such as *linalool*,⁸ *citronellal*,⁹ and *citral*.¹⁰

This preparation is referenced from:

- [Org. Syn. Coll. Vol. 8, 183](#)

References and Notes

1. Department of Synthetic Chemistry, Faculty of Engineering, Shizuoka University, Johoku, Hamamatsu 432 Japan.
2. Chalk, A. J.; Magennis, S. A.; Wertheimer, L. V. in "Fundamental Research In Organometallic Chemistry," Tsutsui, M.; Ishii, Y.; Yaozeng, H., Eds.; Van Nostrand Reinhold: New York, 1982; p. 851.
3. Kofron, W. G.; Baclawski, L. M. *J. Org. Chem.* **1976**, *41*, 1879–1880.
4. Takabe, K.; Katagiri, T.; Tanaka, J. *Tetrahedron Lett.* **1972**, 4009–4012.
5. Takabe, K.; Katagiri, T.; Tanaka, J. *Chem. Lett.* **1977**, 1025–1026.
6. [Stork, G.; Grieco, P. A.; Gregson, M. *Org. Synth., Coll. Vol. VI* **1988**, 638.](#)
7. Sandler, S. R.; Karo, W. "Organic Functional Group Preparations;" Academic Press: New York, 1968, Vol. 1, pp. 324–325.
8. Takabe, K.; Katagiri, T.; Tanaka, J. *Tetrahedron Lett.* **1975**, 3005–3006.
9. Tani, K.; Yamagata, T.; Otsuka, S.; Akutagawa, S.; Kumobayashi, H.; Taketomi, T.; Takaya, H.; Miyashita, A.; Noyori, R. *J. Chem. Soc., Chem. Commun.* **1982**, 600–601.
10. Takabe, K.; Yamada, T.; Katagiri, T. *Chem. Lett.* **1982**, 1987–1988.

Appendix

Chemical Abstracts Nomenclature (Collective Index Number); (Registry Number)

[diethyl ether](#) (60-29-7)

[sodium sulfate](#) (7757-82-6)

[nitrogen](#) (7727-37-9)

[diethylamine](#) (109-89-7)

dimethylamine (124-40-3)

lithium (7439-93-2)

butyllithium (109-72-8)

ISOPRENE (78-79-5)

hexane (110-54-3)

citral (5392-40-5)

calcium hydride (7789-78-8)

Linalool (78-70-6)

neryl chloride

nerol (106-25-2)

citronellal (106-23-0)

N,N-Diethylnerylamine,
2,6-Octadien-1-amine, N,N-diethyl-3,7-dimethyl-, (Z)- (40137-00-6)

dipropylamine (142-84-7)

butyllithium-hexane

N,N-diethyl-2-methyl-2-butenylamine

N,N-diethyl-3-methyl-2-butenylamine

N,N-dimethylnerylamine

N,N-dipropylnerylamine