



A Publication  
of Reliable Methods  
for the Preparation  
of Organic Compounds

## Working with Hazardous Chemicals

The procedures in *Organic Syntheses* are intended for use only by persons with proper training in experimental organic chemistry. All hazardous materials should be handled using the standard procedures for work with chemicals described in references such as "Prudent Practices in the Laboratory" (The National Academies Press, Washington, D.C., 2011; the full text can be accessed free of charge at [http://www.nap.edu/catalog.php?record\\_id=12654](http://www.nap.edu/catalog.php?record_id=12654)). All chemical waste should be disposed of in accordance with local regulations. For general guidelines for the management of chemical waste, see Chapter 8 of Prudent Practices.

In some articles in *Organic Syntheses*, chemical-specific hazards are highlighted in red "Caution Notes" within a procedure. It is important to recognize that the absence of a caution note does not imply that no significant hazards are associated with the chemicals involved in that procedure. Prior to performing a reaction, a thorough risk assessment should be carried out that includes a review of the potential hazards associated with each chemical and experimental operation on the scale that is planned for the procedure. Guidelines for carrying out a risk assessment and for analyzing the hazards associated with chemicals can be found in Chapter 4 of Prudent Practices.


The procedures described in *Organic Syntheses* are provided as published and are conducted at one's own risk. *Organic Syntheses, Inc.*, its Editors, and its Board of Directors do not warrant or guarantee the safety of individuals using these procedures and hereby disclaim any liability for any injuries or damages claimed to have resulted from or related in any way to the procedures herein.

*These paragraphs were added in September 2014. The statements above do not supersede any specific hazard caution notes and safety instructions included in the procedure.*

*Organic Syntheses, Coll. Vol. 8, p.616 (1993); Vol. 66, p.211 (1988).*

## TRISAMMONIUM GERANYL DIPHOSPHATE

[Diphosphoric acid, mono(3,7-dimethyl-2,6-octadienyl) ester (*E*)-, trisammonium salt]

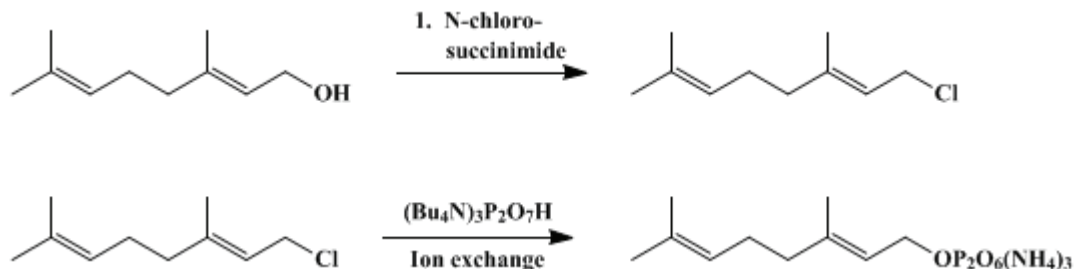

Submitted by Andrew B. Woodside, Zheng Huang, and C. Dale Poulter<sup>1</sup>.

Checked by Pamela Seaton and James D. White.

### 1. Procedure

A. *Geranyl chloride*. To a flame-dried, 100-mL, three-necked, round-bottomed flask equipped with a magnetic stirrer, a low-temperature thermometer, a rubber septum, and a nitrogen inlet adapter is added 1.47 g (11 mmol) of *N*-chlorosuccinimide (Note 1). The powder is dissolved in 45 mL of dry dichloromethane (Note 2), and the resulting solution is cooled to  $-30^{\circ}\text{C}$  with a dry ice–acetonitrile bath. Freshly distilled dimethyl sulfide (0.87 mL, 0.74 g, 12 mmol) is added dropwise by syringe. The mixture is warmed to  $0^{\circ}\text{C}$  with an ice–water bath, maintained at that temperature for 5 min, and cooled to  $-40^{\circ}\text{C}$ . To the resulting milky white suspension is added dropwise by syringe 1.54 g (10 mmol) of geraniol (Note 3) dissolved in 5 mL of dry dichloromethane. The suspension is warmed to  $0^{\circ}\text{C}$  with an ice–water bath and stirred for 2 hr. The ice bath is then removed, and the reaction mixture is allowed to warm to room temperature. Stirring is continued for an additional 15 min. The resulting clear, colorless solution is poured into a 250-mL separatory funnel and washed with 25 mL of saturated sodium chloride. The aqueous layer is washed with two 20-mL portions of pentane. The pentane extracts and an additional 20-mL portion of pentane is added to the methylene chloride extract. The resulting solution is washed twice with 10 mL of saturated sodium chloride and dried over magnesium sulfate. Solid material is removed by vacuum filtration through a fritted-glass funnel, and most of the solvent is removed with a rotary evaporator at aspiratory pressure. The last traces of solvent are removed by pumping at high vacuum (0.2 mm) for 1.5 hr. The resulting pale-yellow oil (1.61 g, 9.3 mmol, 93%) is used directly in the next step (Note 4).

B. *Trisammonium geranyl diphosphate*. To a flame-dried, 100-mL, round-bottomed flask equipped with a magnetic stirrer and a nitrogen inlet adapter is added 9.14 g (9.3 mmol) of tris (tetrabutylammonium) hydrogen pyrophosphate trihydrate (Note 5). The flocculent white solid is dissolved in 20 mL of dry acetonitrile (Note 6). To the resulting milky-white suspension (Note 7) is added 0.83 g (4.8 mmol) of geranyl chloride. The mixture is allowed to stir at room temperature for 2 hr. solvent is then removed with a rotary evaporator using a  $40^{\circ}\text{C}$  water bath. The pale-yellow residue is dissolved in 3 mL of ion-exchange buffer (Note 8), and the resulting clear solution is loaded onto a  $4 \times 15$ -cm column of Dowex AG 50W-X8 (100–200-mesh) cation-exchange resin (ammonium form) (Note 9). The flask is washed twice with 5 mL of buffer and both washes are loaded into the column before elution with 360 mL (two column volumes) of ion-exchange buffer (Note 10). The eluant is collected in a 500-mL freeze-drying flask, frozen as described in Note 5, and lyophilized for 18–24 hr (Note 11) to yield 2.57 g of a white solid. The material is dissolved in 5 mL of 0.05 M ammonium bicarbonate, and the clear solution is transferred to a 50-mL centrifuge tube. Then 20 mL of 1 : 1 (v/v) acetonitrile : isopropyl alcohol is added, and the contents are mixed thoroughly on a vortex mixer, during which time a white precipitate forms. The suspension is cleared by centrifugation for 5 min at 2000 rpm. The supernatant solution is removed with a pipette, and residue is suspended in 5 mL of 0.05 M ammonium

bicarbonate, and the process is repeated. Three additional extractions are performed using 2 mL of 0.05 M ammonium bicarbonate and 8 mL of acetonitrile : isopropyl alcohol. The combined supernatant solutions (approximately 80 mL) are concentrated to approximately 5 mL with a rotary evaporator at 40°C (Note 12).

One-half of the concentrated extract, dissolved in an equal volume of chromatography buffer (Note 13), is loaded onto a 5.5 × 18-cm cellulose flash column<sup>2</sup> (Note 14). The flask is rinsed with three 5-mL portions of chromatography buffer and each is loaded onto the column. The column is then eluted with 900 mL of chromatography buffer. After a 50-mL forerun, twenty-eight 30-mL fractions are collected, and every second fraction is analyzed by thin-layer chromatography (Note 15). Fractions containing trisammonium geranyl diphosphate (typically 12–23) are pooled and concentrated to approximately 120 mL with a rotary evaporator at 40°C. The concentrate is transferred to a 250-mL freeze-drying flask and lyophilized for 18–24 hr as previously described in (Note 5). The resulting flocculent white solid is collected and stored at –78°C. The cellulose chromatography is repeated to yield a total of 1.51–1.55 g (85–87%) of trisammonium geranyl diphosphate from geraniol (Note 16).

## 2. Notes

1. *N*-Chlorosuccinimide (from the Aldrich Chemical Company, Inc.) is recrystallized from benzene (*Caution: Carcinogenic*).
2. Methylene chloride is distilled from phosphorus pentoxide immediately before use.
3. Geraniol (from the Aldrich Chemical Company, Inc.) is distilled before use, bp 90–92°C at 3 mm.
4. The IR, <sup>1</sup>H NMR, and <sup>13</sup>C NMR spectra of this material are identical to those for distilled geranyl chloride (bp 49–51°C at 0.2 mm). Distillation on a small scale significantly reduces the yield, and there is no improvement in the yield of the phosphorylation reaction using distilled material. A synthesis of geranyl chloride was reported earlier in this series.<sup>3</sup> We find, however, that the procedure of Corey et al.<sup>4</sup> is more convenient.
5. Disodium dihydrogen pyrophosphate (3.13 g, 14 mmol) (from Sigma Chemical Co.) is dissolved in 25 mL of deionized water containing 1 mL of concentrated ammonium hydroxide. The resulting clear solution is loaded onto a 2 × 30-cm column of Dowex AG 50W-X8 cation exchange resin (100–200 mesh, H<sup>+</sup>) and eluted with deionized water. The first 150 mL of eluant is collected in a 250-mL freeze-drying flask. A magnetic stirring bar is added, and the solution is titrated to pH 7.3 by slow addition of tetrabutylammonium hydroxide (Aldrich Chemical Company, Inc.). The stirring bar is removed, and the flask is placed in a dry ice–propanol bath. The flask is spun slowly in a manner to uniformly freeze its contents to the walls. Water is removed by lyophilization for 24 hr. The resulting flocculant white solid (10.6 g, 83%) contains three to four waters of hydration and is used without further purification. The material is extremely hygroscopic and can be stored in a desiccator over phosphorus pentoxide.
6. Reagent-grade acetonitrile is distilled from phosphorus pentoxide immediately before use.
7. Clear solutions can be obtained by filtration. Residual water can be removed from freeze-dried salt by repeated evaporation of dry acetonitrile (rotary evaporator). This material is not noticeably more effective than the salt obtained after lyophilization.
8. Ammonium bicarbonate (2.0 g) is dissolved in 1 L of 2% (v/v) isopropyl alcohol/water. The resulting solution is 25 mM in ammonium bicarbonate.
9. The ammonium form of the resin is generated by placing 188 mL of Dowex AG 50W-X8 (100–200 mesh, H<sup>+</sup> form) in an 1-L fritted-glass funnel and washing the material with four 200-mL portions of concentrated ammonium hydroxide. The resin is washed with 200-mL portions of deionized water until the pH of the filtrate drops to pH 7, then with two 200-mL portions of ion-exchange buffer. The washed resin is suspended in 200 mL of buffer and slurry-packed into the column.
10. Dowex AG 50W-X8 (100–200 mesh) from BioRad has a capacity of 1.7 meq per milliliter of resin bed. This represents approximately 10-fold excess of exchangeable ions in the resin over material loaded onto the column. However, the tetrabutylammonium cation has a lower affinity for the resin than the ammonium cation. To optimize the exchange, it is important to maintain a low concentration of ammonium ion in the exchange buffer, to elute the material slowly (less than 9 mL/min on the 4 × 15-cm column), and to elute with only two column volumes of exchange buffer; otherwise previously exchanged tetrabutylammonium cation will begin to elute from the column. Incomplete exchange dramatically reduces the efficiency of the subsequent purification on cellulose. The efficiency of the exchange can be determined by <sup>1</sup>H NMR.

11. **Trisammonium geranyl diphosphate** will decompose if left under vacuum for extended periods. It is important to remove the sample from the freeze drier within a few hours after water has been removed.
12. This material is stored at  $-20^{\circ}\text{C}$  until chromatography on cellulose.
13. Chromatography buffer is prepared by dissolving 4.0 g of **ammonium bicarbonate** in 250 mL of deionized water and adding 500 mL of **isopropyl alcohol** and 250 mL of **acetonitrile**. The resulting solution is approximately 50 mM in **ammonium bicarbonate**.
14. Whatman CF11 fibrous cellulose powder is prepared for chromatography by the following procedure. Cellulose powder (1 L, dry volume) is mixed with 700 mL of deionized water in a 2-L beaker by vigorous stirring with a glass rod. The suspension is allowed to stand for 30 min, and the water is removed by decantation. The same procedure is followed as the cellulose is washed in succession with two 700-mL portions of 0.1 N **hydrochloric acid**, two 700-mL portions of deionized water, two 700-mL portions of 0.1 N **sodium hydroxide**, two 700-mL portions of deionized water, and two 700-mL portions of 1 : 1 (v/v) **isopropyl alcohol** : water. The material is stored at  $4^{\circ}\text{C}$  in 1 : 1 **isopropyl alcohol** : water until used. The column is slurry-packed in 1 : 1 **isopropyl alcohol** : water and washed with 1.3 L (approximately three column volumes) of **acetonitrile**. The column is then washed with 1.3 L of 1 : 1 **isopropyl alcohol** : water and equilibrated with 1.3 L of chromatography buffer (Note 13).
15. E. Merck cellulose thin-layer chromatography plates (available from American Scientific Products) are developed with chromatography buffer (Note 13) and visualized with **sulfosalicylic acid/ferric chloride** spray.<sup>5</sup> The system consists of a solution of 1.0 g of **sulfosalicylic acid** (from Aldrich Chemical Co., Inc.) dissolved in 100 mL of 3:2 (v/v) **ethanol** : water and a solution of 0.20 g of **ferric chloride** in 100 mL of 4:1 (v/v) **ethanol** : water. Plates are first sprayed with the **sulfosalicylic acid** solution (thoroughly wetted but not dripping) and allowed to air-dry. The **ferric chloride** solution is lightly sprayed onto the plates. Phosphate-containing compounds appear as white spots on a pink background. A second light spraying with **ferric chloride** may be necessary to make the spots pronounced. It is *important* to prepare both sprays freshly. Their shelf life is only about 6 hr. Under these conditions the **trisammonium geranyl diphosphate** has an  $R_f$  of 0.35. Residual tetrabutylammonium salt moves with the solvent front, and ammonium inorganic pyrophosphate remains at the origin.
16. This material migrates as a single spot on the cellulose thin-layer system and has no extraneous peaks in the  $^1\text{H}$ ,  $^{13}\text{C}$ , and  $^{31}\text{P}$  NMR spectra. The IR and NMR spectral properties of **trisammonium geranyl diphosphate** are as follows: IR (KBr)  $\text{cm}^{-1}$ : 3100–3500 (br), 2800, 2990 (br), 1650, 1450, 1400, 1200, 1120, 1080, 1015, 900;  $^1\text{H}$  NMR (300 MHz,  $\text{D}_2\text{O}/\text{ND}_4\text{OD}$ )  $\delta$ : 1.62 (s, 3 H, methyl), 1.68 (s, 3 H, methyl) 1.72 (s, 3 H, methyl), 2.11 (m, 4 H,  $\text{CH}_2$  at  $\text{C}_4$  and  $\text{C}_5$ ), 4.47 (t, 2 H,  $J_{\text{1H,1H}} = 6.5$ ,  $J_{\text{1H,31P}} = 6.5$ ,  $\text{CH}_2$  at  $\text{C}_1$ ), 5.22 (broad, 1 H,  $J_{\text{1H,1H}} = 6.5$ , H at  $\text{C}_6$ ), and 5.47 (t, 1 H,  $J_{\text{1H,1H}} = 6.5$ , H at  $\text{C}_2$ );  $^{13}\text{C}$  NMR (75 MHz,  $\text{D}_2\text{O}/\text{ND}_4\text{OD}$ ,  $^1\text{H}$  decoupled): 18.28 ( $\text{CH}_3$ ), 19.69 ( $\text{CH}_3$ ), 27.55 ( $\text{CH}_3$ ), 28.60 ( $\text{CH}_2$ ), 41.55 ( $\text{CH}_2$ ), 65.42 ( $\text{CH}_2$ , d,  $J_{\text{13C,31P}} = 4.0$ ), 122.85 (CH, d,  $J_{\text{13C,31P}} = 7.5$ ), 127.10 (CH), 136.68 (C), 145.76 (C);  $^{31}\text{P}$  NMR (32 MHz,  $\text{D}_2\text{O}/\text{ND}_4\text{OD}$ ,  $^1\text{H}$  decoupled):  $-11.23$  (d, 1 P  $J_{\text{31P,31P}} = 20$ ) and  $-9.10$  (d, 1 P,  $\text{P}_2$ ).

### 3. Discussion

Previous methods for the preparation of salts of **geranyl diphosphate** and other allylic isoprenoid diphosphates are based on condensation between the alcohol and inorganic phosphate by **trichloroacetonitrile** as originally reported by Cramer<sup>6</sup> and modified by Popjak.<sup>7</sup> The reaction generates a complex mixture of organic and inorganic polyphosphates that must be separated by chromatography. The desired diphosphate ester has been prepared on small scales in yields of up to 30%,<sup>8</sup> but in our experience, the yields of pure material obtained by this procedure are usually less than 10%.

The direct displacement reaction can be used to prepare many of the common diphosphate esters in the **isoprene** biosynthetic pathway,<sup>9</sup> including **isopentenyl diphosphate**.<sup>10,11</sup> The yields are typically 60–90% from the alcohol, and the absence of phosphate polymer found in the Cramer procedure simplifies the purification step. We have also used the displacement procedure to prepare radio-labeled material for biosynthetic studies.<sup>12</sup>

---

## References and Notes

1. Department of Chemistry, University of Utah, Salt Lake City, UT 84112.

2. Still, W. C.; Kahn, M.; Mitra, A. *J. Org. Chem.* **1978**, *48*, 2923.
  3. Calzada, J. G.; Hooz, J. *Org. Synth., Coll. Vol. VI* **1988**, 634.
  4. Corey, E. J.; Kim, C. U.; Takeda, M. *Tetrahedron Let.* **1972**, 4339.
  5. "CRC Handbook of Chromatography," Gunter, Z.; Sherma, J. Eds.; CRC Press: Cleveland, OH, 1972, Vol. 2, p. 143.
  6. Cramer, F.; Rittersdorf, W.; Bohm, M. *Liebigs Ann. Chem.* **1962**, *654*, 180.
  7. Popjak, G.; Cornforth, J. W.; Cornforth, R. H.; Rhyage, R.; Goodman, D. S. *J. Biol. Chem.* **1962**, *237*, 56.
  8. Banthorpe, D. V.; Christou, P. N.; Pink, C. R.; Watson, D. G. *Phytochemistry* **1983**, *22*, 2465.
  9. Dixit, V. M.; Laskovics, F. M.; Noall, W. I.; Poulter, C. D. *J. Org. Chem.* **1981**, *46*, 1967.
  10. Davisson, V. J.; Woodside, A. B.; Poulter, C. D. *Methods Enzymol.* **1985**, *110*, 130–144.
  11. Davisson, V. J.; Woodside, A. B.; Neal, T. R.; Stremmer, K. E.; Muehlbacher, M.; Poulter, C. D. *J. Org. Chem.* **1986**, *51*, 4768–4779.
  12. Davisson, V. J.; Zabriskie, T. M.; Poulter, C. D. *Bioorg. Chem.* **1986**, *14*, 46–54.
- 

**Appendix**  
**Chemical Abstracts Nomenclature (Collective Index Number);**  
**(Registry Number)**

Diphosphoric acid, mono(3,7-dimethyl-2,6-octadienyl) ester (E)-, trisammonium salt

ethanol (64-17-5)

hydrochloric acid (7647-01-0)

Benzene (71-43-2)

acetonitrile (75-05-8)

sodium hydroxide (1310-73-2)

sodium chloride (7647-14-5)

isopropyl alcohol (67-63-0)

ferric chloride (7705-08-0)

ammonium hydroxide (1336-21-6)

sulfosalicylic acid

N-chlorosuccinimide (128-09-6)

Pentane (109-66-0)

methylene chloride,  
dichloromethane (75-09-2)

trichloroacetonitrile (545-06-2)

magnesium sulfate (7487-88-9)

dimethyl sulfide (75-18-3)

ISOPRENE (78-79-5)

geraniol (106-24-1)

tetrabutylammonium hydroxide (2052-49-5)

Geranyl chloride (5389-87-7)

phosphorus pentoxide (1314-56-3)

Trisammonium geranyl diphosphate

hydrogen pyrophosphate trihydrate

ammonium bicarbonate

Disodium dihydrogen pyrophosphate

geranyl diphosphate

isopentenyl diphosphate